

Daniel

TG-142

Linac Calibration & QA Laboratory 2025

Exercise 1: Optical distance indicator vs front pointer

Questions & Answers:

1. Describe how we can test the optical distance indicator and any assumptions inherent in this test

ODI

Front pointer

We use a metallic calibrated rod to match the actual wanted distances with the ODI

2. What did you observe when you performed this test?

The front pointer needs a movement
there's different size front pointers 90cm, 100cm, 110cm
adjust to the surface
±1mm in ODI
turn off lights, look ODI, is it at the desired d? ✓

- If not...
3. What can the clinical physicist do if he or she detects an error in the optical distance indicator?

check tolerance
based on
TG142
if off

call the engineer

100SSD being the
most important to
check

Exercise 2: Laser positioning and isocentricity

Questions & Answers:

4. What is the tolerance for laser localization?

1mm

5. Why is laser localization important? (Think about what could be impacted clinically if your lasers are off)


The radiologist uses the lasers as guide and on daily QA we trust these, so if they are off we get off sets on dose profiles.

6. What can the clinical physicist do if he or she detects an error in laser positioning?

call the engineer

Again only if it is truly off tolerance

7. What is the tolerance for isocenter coincidence?

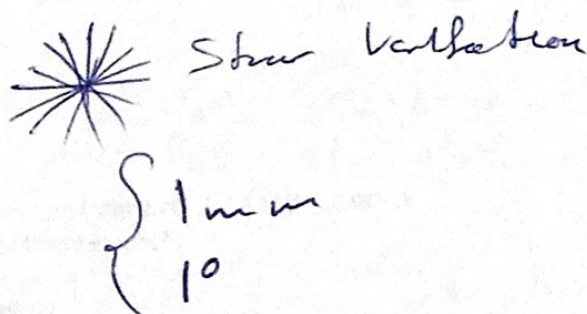

Exercise 3: Radiation light field coincidence


Questions & Answers:

8. What is the tolerance for radiation field light coincidence?

3% Daily
2% Monthly
1% Annual

9. Why is coincidence of the radiation field and the light field important?

This ensures the light field is accurately representing the radiation field during treatment

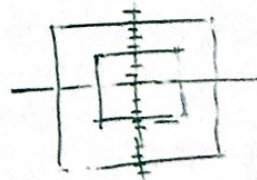

- The light field is used during setup of both phantoms and patients

10. What can the clinical physicist do if there is an error in the radiation field light coincidence?

call the engineer

OPI
Shield light } inside gully head
Lasers in the well.

Exercise 4: Linac output calibration (TG51)

Questions & Answers:

11. Write down your values for P_{TP} , P_{pol} , P_{ion} . What are typical ranges of values for P_{TP} , P_{pol} , and P_{ion} ?

$P_{TP} = 1.0786$? $32^{\circ}C$
 $P_{pol} = 1.002$ ✓
 $P_{ion} = 1.0015$

note \rightarrow Electrometer readings
correction factors

around 7 If $\sim 22^{\circ}C$ & $101.33 kPa$
 P_{TP}

P_{pol} 0.996 - 1.004

P_{ion} 1.0 - 1.05

12. What kinds of errors can contribute to the TG51 result deviating from 1 cGy/MU at d_{max} ? What can the clinical physicist do when he or she detects each of these errors?

Setup Errors \rightarrow Phantom Chamber

OPI
light field / cross hair

Verify source of error, consider recent Daily QA

13. Why are most of our TG51 measurements obtained at 10cm instead of, say, d_{max} ?

10 cm is reproducible, reliable

d_{max} you have to estimate
and changes or errors in this,

~~scale up more~~
have a bigger impact


Additional questions:

14. How do each of the parameters tested relate to one another (e.g., how does laser localization impact TG51?), and what order would you perform the tests from this lab in the clinic?

Errors will cascade into the calculations.

If ODS, field size, or any other "mechanical" reference is off this means all calculations are off. \rightarrow test all dependent measurements sequentially at the end

15. Describe how an ionization chamber works

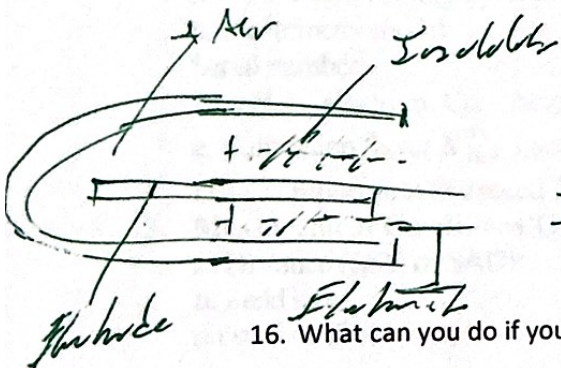

A voltage is applied to polarize ions and capture them and obtain their charges or current

Ionizing radiation produces these ions on the air inside. the Electric field between V pull them out for measurement

16. What can you do if your chamber is not listed in TG51 or the addendum to obtain a k_Q value?

You make sure what you have matches physical properties of one listed. Institutions should get you the right chamber.

17. What is the purpose of k_Q ?

Calibration factor for reference dosimetry
It takes values from ^{60}Co source to our setup

18. How do you know if your chamber can be used for "reference dosimetry"?

It is listed in TG51 and complies with requirements in table III

TG-51 Worksheet A: Photon Beams

1. Site data

Institution:

Physicist:

Date:

Accel or ^{60}Co Mfr.:

Model & serial number:

Nominal photon energy/beam identifier:

2. Instrumentation

a. Chamber model:

Serial number:

Cavity inner radius (r_{cav} , Table III):

Waterproof:

If no, is waterproofing ≤ 1 mm PMMA or thin latex?:

b. Electrometer model:

Serial number:

i. P_{elec} , electrom. Corr factor (Sec. VII.B):

c. Calibration factor $N_{D,w}^{60\text{Co}}$ (Sec.V):

Date of report (not to exceed 2 years):

3. Measurement Conditions ($10 \times 10 \text{ cm}^2$, point of measurement at 10 cm depth (water equivalent))

a. Distance (SSD or SAD):

b. Field size:

on surface (SSD setup)

at detector (SAD setup):

c. Number of monitor units:

4. Beam Quality (Sec. VIII.B – not needed for ^{60}Co)

If energy < 10 MV, use no lead foil.

Measure %dd(10) [% depth-dose at 10 cm depth for curve shifted upstream by $0.6r_{\text{cav}}$]

Field size $10 \times 10 \text{ cm}^2$ on surface, SSD = 100 cm: yes

a. %dd(10)_x = % dd(10)

☒ If energy ≥ 10 MV

Distance of 1 mm lead foil phantom surface

Measure % dd(10)_{Pb} [% depth-dose at 10 cm depth for curve shifted upstream by $0.6 r_{\text{cav}}$]

Field size $10 \times 10 \text{ cm}^2$ on surface, SSD = 100 cm: yes

%dd(10)_{Pb} (includes e^- contamination):

50 cm: %dd(10)_x = $[0.8905 + 0.00150\%dd(10)_{\text{Pb}}]\%dd(10)_{\text{Pb}}$ [%dd(10)_{Pb} $\geq 73\%$] Eq.(13)

30 cm: %dd(10)_x = $[0.8116 + 0.00264\%dd(10)_{\text{Pb}}]\%dd(10)_{\text{Pb}}$ [%dd(10)_{Pb} $\geq 71\%$] Eq. (14)

If %dd(10)_{Pb} $< 71\%$ (30 cm) or 73% (50 cm): %dd(10)_x = %dd(10)_{Pb}

b. %dd(10)_x (for open beam):

Has lead foil been removed?

☒ Interim alternative for energy > 10 MV & with ≥ 45 cm clearance: using no lead foil

Measure %dd(10) [% depth-dose at 10 cm depth for curve shifted upstream by $0.6 r_{\text{cav}}$]

%dd(10): No lead

%dd(10)_x = $1.267 (\%dd(10) - 20.0)$ [for $75\% < \%dd(10) \leq 89\%$]

c. %dd(10)_x =

TG-51 Worksheet A: Photon Beams (cont'd)

MU = 200

5 Determination of k_Q (Sec. IX. B)

Chamber model used to get k_Q :

a. %dd(10)_x (from 4, above):

b. k_Q [Table I or Fig. 4]:

$$\frac{PTW - N3103}{66.2\%} = 0.992$$

6 Temperature/pressure Correction (Sec. VII.C)

a. Temperature:

b. Pressure:

c. P_{TP} :

$$\frac{32.2}{28.7 \text{ in/Hg}} \times \frac{101.33}{760} = 1.0786$$

$$97.19 \text{ kPa} = 728.98 \text{ mmHg}$$

$$Eq.(10) = \left(\frac{273.2 + 6a}{295.2} \right) \left(\frac{101.33}{6b} \right)$$

7 Polarity correction (Sec. VII. A.)

M_{raw}^+ :

M_{raw}^- :

a. M_{raw} (for polarity of calibration):

b. P_{pol} :

$$\frac{12.65}{12.7} = 1.002$$

$$Eq.(9) = \left| \frac{M_{raw}^+ - M_{raw}^-}{2M_{raw}} \right|$$

8 P_{ion} measurements (Sec. VII. D. 2)

Operating voltage = V_H :

Lower voltage V_L :

M_{raw}^H :

M_{raw}^L :

^{60}Co treated as general recombination

a. $P_{ion}(V_H)$ (Eq.(11)):

$$\frac{300}{150} = 2$$

$$\frac{12.65}{12.63} = 1.0005$$

$$\left[\left(1 - \left(\frac{V_H}{V_L} \right)^2 \right) / \left(\frac{M_{raw}^H}{M_{raw}^L} - \left(\frac{V_H}{V_L} \right)^2 \right) \right]$$

Pulsed/swept beams

b. $P_{ion}(V_H)$ (Eq.(12))

$$1.0015 \left[\left(1 - \frac{V_H}{V_L} \right) / \left(\frac{M_{raw}^H}{M_{raw}^L} - \frac{V_H}{V_L} \right) \right]$$

If $P_{ion} > 1.05$, another ion chamber should be used.

9. Corrected ion. ch. rdg. M (Sec. VII) at 10 cm depth, water equivalent

$$M = P_{ion} P_{TP} P_{elec} P_{Pol} M_{raw} = [8(a \cdot \text{or} \cdot b) \cdot 6c \cdot 2bi \cdot 7b \cdot 7a]$$

Fully corrected M (Eq.(8)):

10. Dose to water at 10 cm depth:

a. Dose to water at 10 cm depth =

b. Dose/ MU(or min ^{60}Co) at 10 cm depth

$$D_w^Q = Mk_Q N_{D,w}^{60Co} = [9.5b \cdot 2c] \text{ Eq. (3)}$$

$$1.371822 \text{ Gy}$$

$$6.859112 \text{ Gy/MU} [10a/3c]$$

11. Dose to water/MU(or min, ^{60}Co) at d_{max} (if relevant locally)

a. Clinical %dd(10) for SSD setup / 100.:
or clinical TMR(10, 10 x 10) for SAD setup:

b. Dose / MU(or min, ^{60}Co) at d_{max} :

$$66.2\%$$

$$1.036119 \text{ cGy (water)/MU}$$